


Subperiosteal implants

A multidisciplinary team plans, designs and manufactures the customized structure for the subperiosteal implant treatment.

Its individualized design maximizes stability and facilitates bone regeneration providing immediate esthetics and functional results.

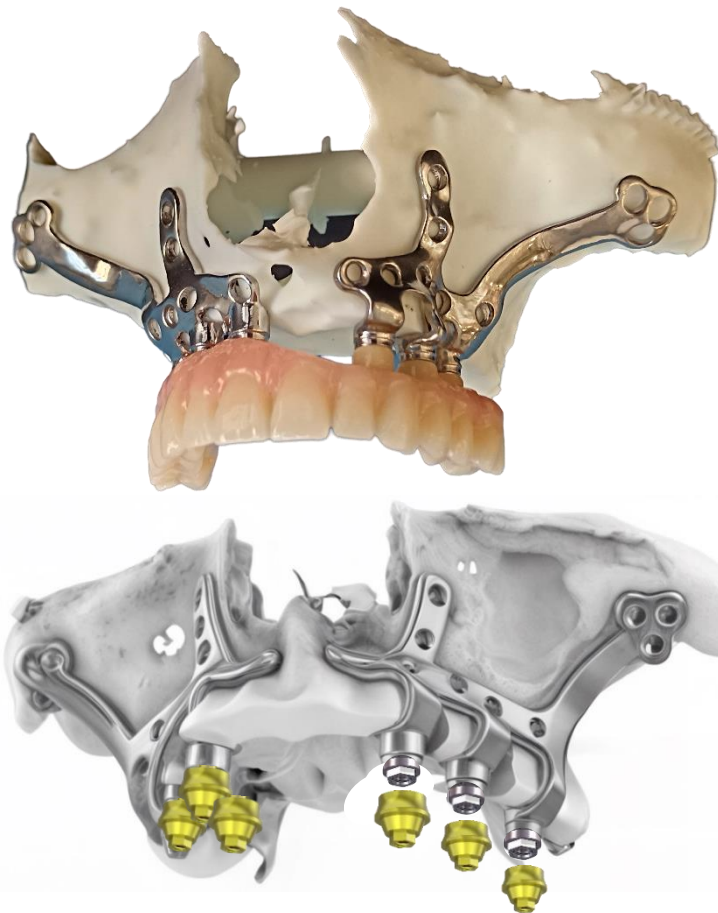


Advantages


- Multi-Unit removable abutments of different heights. Nitrided in yellow (PVD).
 - Helical CT or Cone Beam thanks to post-processing with artificial intelligence.
 - Titanium framework (Ti-6Al-4V), direct laser sintered, SLM.
 - Framework thickness 0.8 mm.
 - Full arch treatments with independent frameworks per quadrant and positioning guide.
- Extensive international experience with more than 250 subperiosteal implant treatments performed annually.
 - Specific planning software for this type of treatment (Rhino Medical).
 - Delivery of all necessary instruments for surgery and assistance of a product specialist.

With the collaboration of

KUNE
IMPLANTS

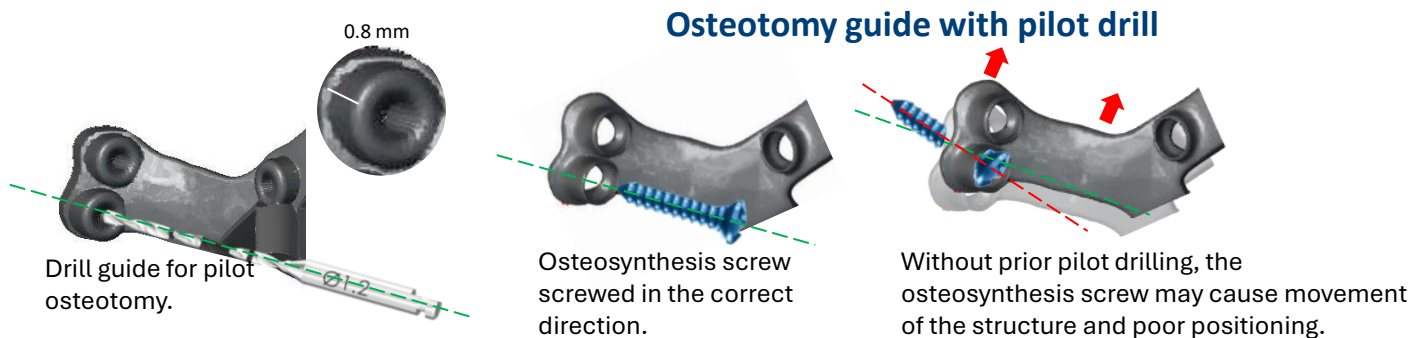


Osteotomy guide with osteosynthesis screw template

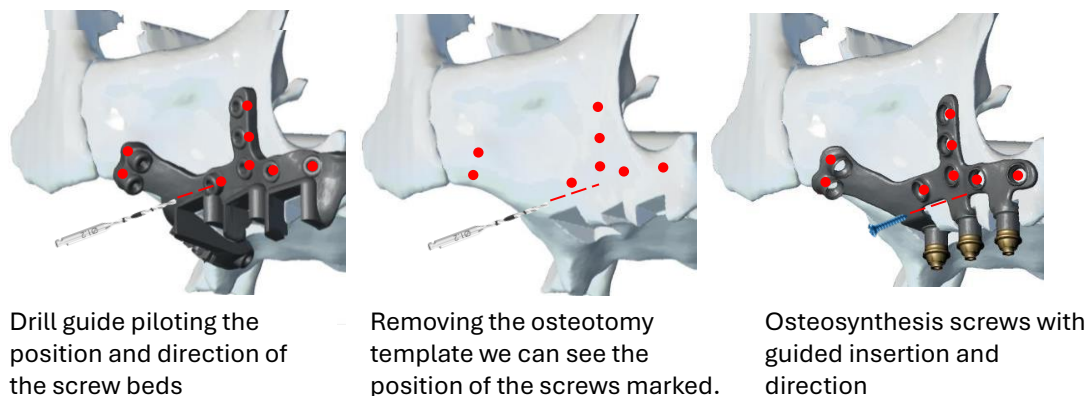

Osteotomy drill Ø 2mm




Design of the osteotomy guide drawers

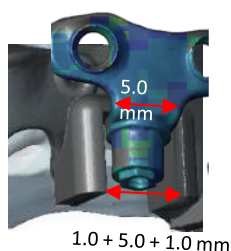


We added 1mm extra per side in the guide "drawers" for a correct milling, which allows a better seating of the arms of the structures.

Macro osteotomy guide

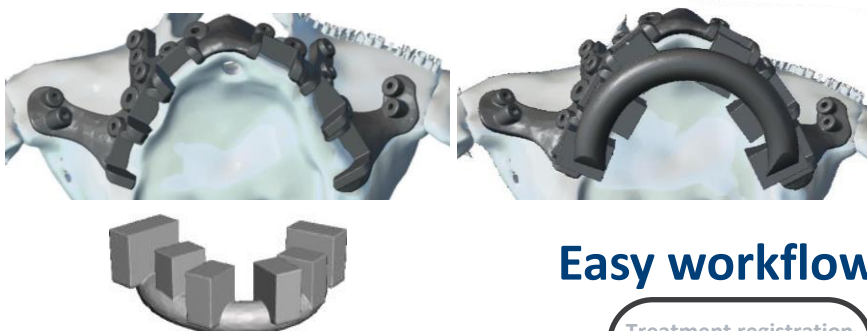

Enlarged vestibular step

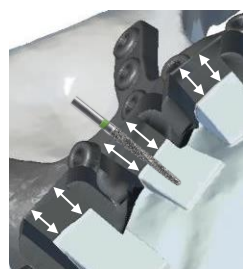

Wider vestibular step that provides a support point for the osteotomy drill and aids in its correct direction



The macro structure fits with the osteotomy guide to verify correct osteotomies of the drawers.

Pack subperiosteal implants

- Conversion of double CT Cone Beam (clinic) or helical (hospital).
- Design and planning of the subperiosteal frameworks by engineers specialized in craniofacial and maxillofacial medicine.
- Sending the planning by e-mail.
- Review and modification in remote control with the engineers.
- Creation of the surgical report.
- Fabrication and shipment of the subperiosteal implants, cutting guide and positioning guide if necessary.
- Delivery to the clinic of all the surgical material, osteosynthesis screws, rescue screw, screwdriver tips, handle and osteosynthesis drills.



Easy workflow

